

UNIQUE COUPLING BETWEEN μ DSC AND VOLUMETRY

CO₂ ADSORPTION ON BASOLITE C300 STUDIED BY COUPLED VOLUMETRY - MDS

Metal-organic frameworks (MOFs) offer an advantage over other classical porous materials (activated carbon, zeolites) because their properties can be tailored for specific applications. Today these structures are envisaged for multiple applications in gas adsorption, like hydrogen storage, selective adsorption of CO₂ against CH₄, CO₂ against H₂. The Cu-BTC (fig 2.) is chosen as a example of Metal Organic Framework that is available commercially to demonstrate adsorption studied by Sievert’s technique and heat flow calorimetry.

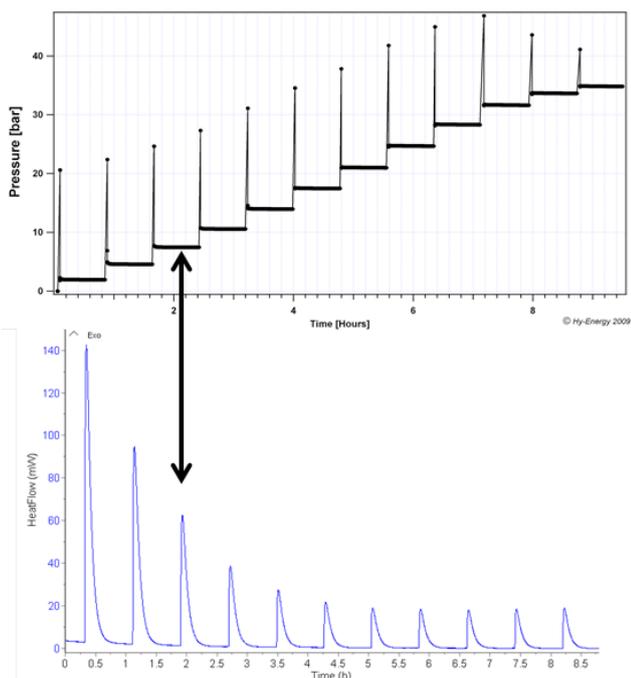


Figure 1. Up: Pressure versus time during a CO₂ adsorption experiment. At the beginning of each aliquot, the step is the measurement of the pressure in the manifold. Down: Calorimetric Heat Flow signal versus time. Each exothermic peak corresponds to the adsorption of a gas dose.

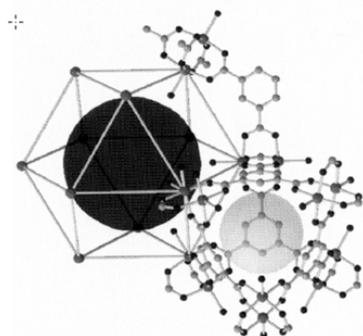


Figure 2. Basolite C300 Cu-BTC

EXPERIMENTAL

CO₂ adsorption into Basolite C300 was measured at -20°C, 30°C and 50°C using a PCTPro-E&E Sievert’s apparatus coupled with a μ DSC7 evo. Gas density temperature correction was performed by measuring the apparent free gas volume at these temperatures.

Understanding the thermodynamics of the adsorption is essential for the practical application and among all the enthalpy of adsorption (or desorption) it is a key parameter.

There are two ways to determine this enthalpy:

- The first one is an indirect method, where the enthalpy is derived from adsorption isotherms at different temperatures.
- The second one is a direct method, where the enthalpy is measured via a calorimetric technique. The biggest disadvantage of this technique is that it gives a result per mole of solid sample and not per mole of gas.

The combination of manometric technique (to quantify the amount of hydrogen absorbed/ released) and calorimetry overcomes this issue and the direct measurement of enthalpy of formation per mole of gas is presented here (figure 1.).

KEY WORDS
ADSORPTION
CO₂
MOFS
VOLUMETRY
CALORIMETRY



μ DSC7 evo
 -45°C to 120°C



sales@setaram.com
www.setaram.com

UNIQUE COUPLING BETWEEN μ DSC AND VOLUMETRY

RESULTS AND CONCLUSIONS

Coupling allows the access to a full characterization of the adsorption: adsorption isotherms (fig 3.), integral heats of adsorption (fig 4.), isotheric heat of adsorption (fig 5.), differential heat of adsorption as a function of the adsorbed amounts (fig 6.).

The isosteric heat of adsorption is obtained by indirect methods by using the isotherms at different temperature. It gives the heat of adsorption at a specific coverage. It is obtained according to:

$$q_{st} = -RT^2 \left(\frac{\partial \ln P}{\partial T} \right)_\theta$$

The differential heat of adsorption is the energy release during the adsorption for a infinitely small quantity of adsorbent. It depends on the adsorbed quantity of the gas. Practically this value is very important. It is this value that is used for the calculation of the sorption installations.

$$q_{diff} = \frac{\partial q_{int}}{\partial n^a}$$

The PCTPro and μ DSC7 evo combination is an ideal technique for the detailed characterization of high surface area materials used in gas sorption. The ease of use and the temperature and pressure ranges are perfect for this type of material application and in particular in the case of CO₂ sorption studies.

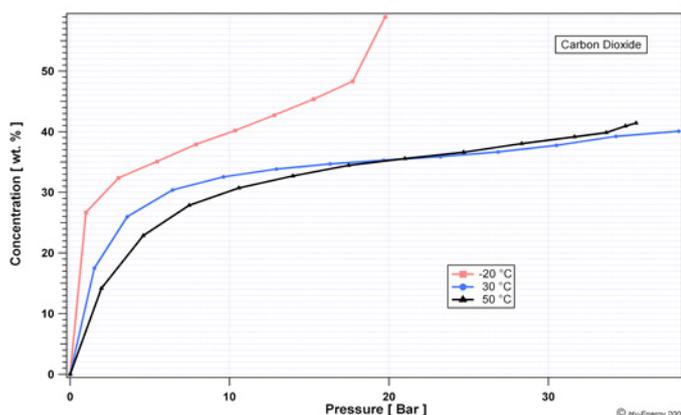


Figure 3. CO₂ adsorption isotherms at -20, 30, 50 °C

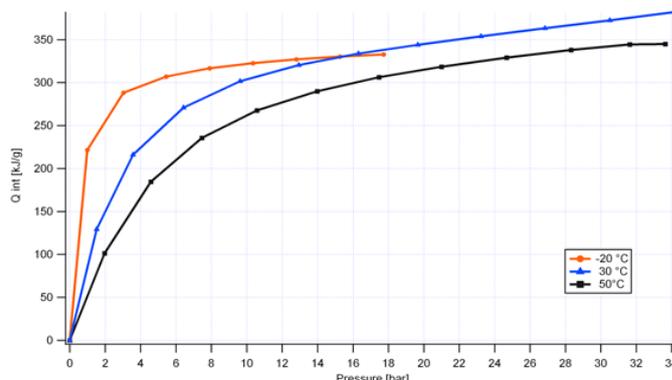


Figure 4. Total heat isotherms at -20, 30, 50 °C. It is the sum of all heat of adsorption versus pressure 30, 50 °C

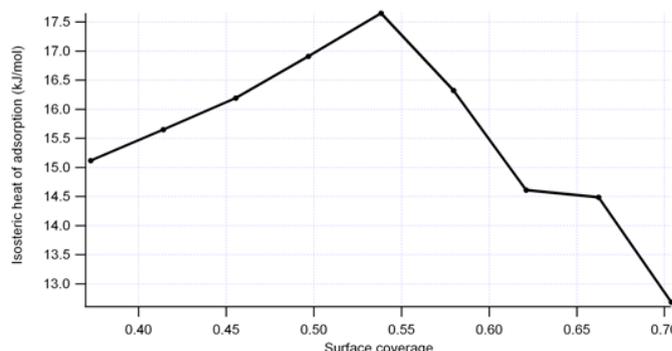


Figure 5. Isotheric heat of adsorption versus surface coverage

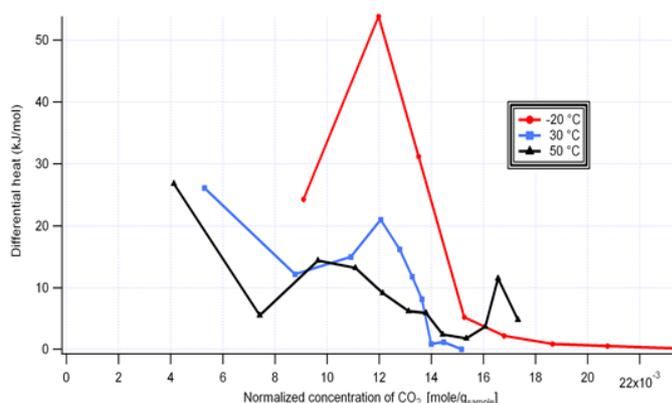


Figure 6. Differential heat of adsorption versus the normalized CO₂ concentration

PCTPro
-260°C to 500°C

KEY WORDS
ADSORPTION
 CO₂
 MOFS
VOLUMETRY
CALORIMETRY

SETARAM
Instrumentation
K&E Technologies

sales@setaram.com
www.setaram.com